

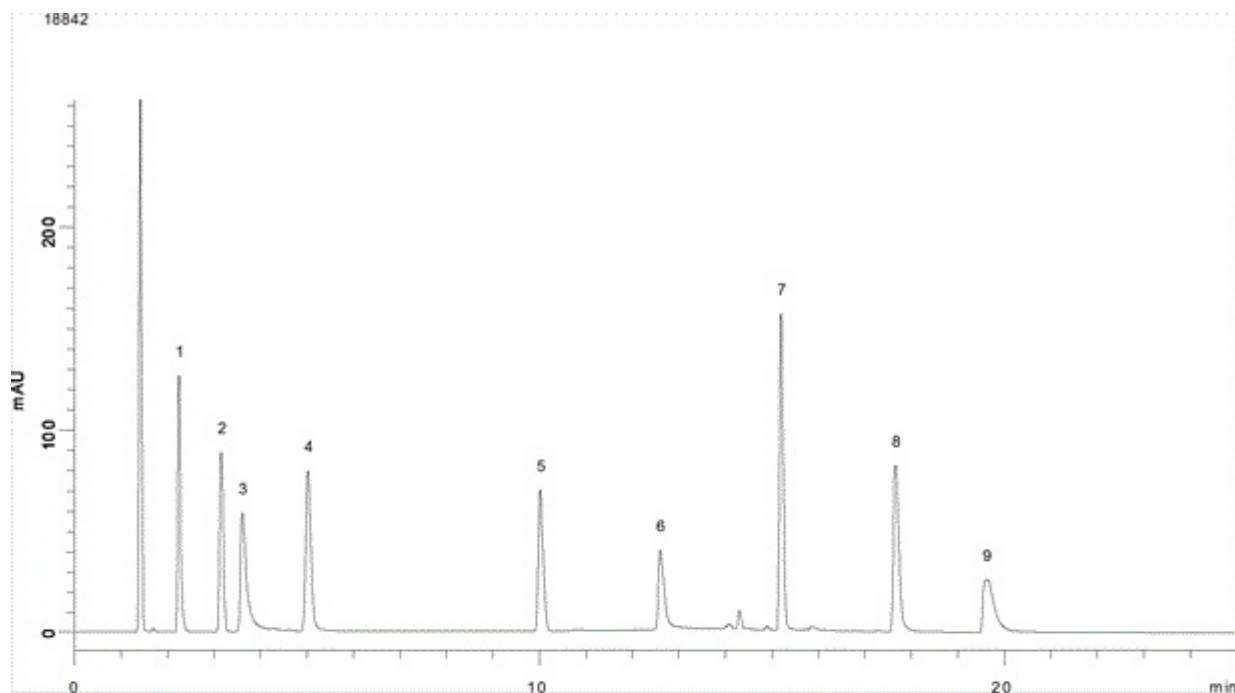


Figure 2: Water-soluble vitamins on Kinetex 2.6 um HILIC 150 x 4.6 mm.

A: Acetonitrile

B: 100 mM Ammonium formate pH 3.2

C: Water

Time (min)	Pct A	Pct B	Pct C
0	90	10	0
5	90	10	0
20	50	10	40

1 mL/min, ambient temperature, UV-Vis 260 nm

1: Nicotinamide

2: PABA

3: Pyridoxine

4: Riboflavin

5: Nicotinic acid

6: Ascorbic acid

7: Folic acid

8: Vitamin B12

9: Thiamine

Caption:

Description:

Dimensions: 613 x 727

aperture: 0

credit:

camera:

caption:

created_timestamp: 0

copyright:

focal_length: 0

iso: 0

shutter_speed: 0

title:

orientation: 0

keywords: Array