


Figure 4. Phenomenex Luna® Omega 1.6µm C18 50 x 2.1mm

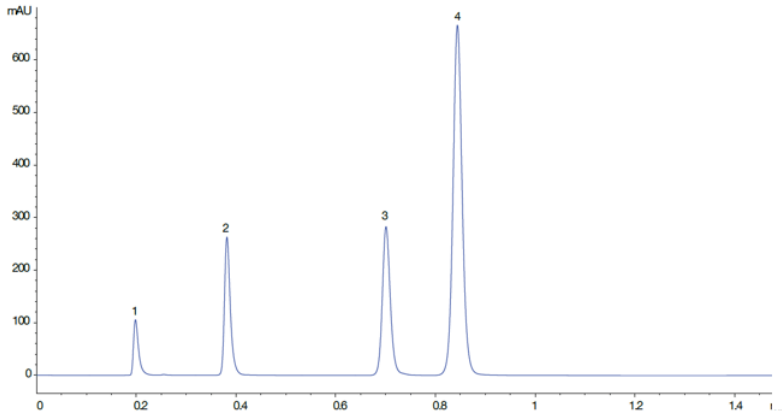

UHPLC Conditions for all Columns:

Column: Luna Omega 1.6 µm C18
Agilent ZORBAX Eclipse Plus 1.8 µm C18 RRHD
Waters ACQUITY UPLC 1.7µm BEH C18
Dimensions: 50 x 2.1 mm
SecurityLINK Part No.: AJ1-1421, AJ1-1441
Mobile Phase: Water/Acetonitrile (35:65)
Flow Rate: 0.5 mL/min
Temperature: 25°C
Injection Volume: 0.2 µL
Sample:

1. Uracil
2. Acetophenone
3. Toluene
4. Naphthalene

Caption:

Description:

Dimensions: 1133 x 432

aperture: 0

credit:

camera:

caption:

created_timestamp: 0

copyright:

focal_length: 0

iso: 0

shutter_speed: 0

title:

orientation: 0

keywords: Array